

Fast, Accurate, Easy to Use

Gluco Spark[®]

Blood Glucose Measuring System



Features:

- ▶ 0.5 μ L Blood Sample Volume.
- ▶ 5 Seconds Fast Test Time.
- ▶ No Coding.
- ▶ Large LCD Display.
- ▶ Strip Ejection.

Life Time
WARRANTY

CE IVD


ISO
15197 - 2013

ISO
9001 - 2015

ISO
13485 - 2016

SENSA  CORE[®]

Strip Advantages

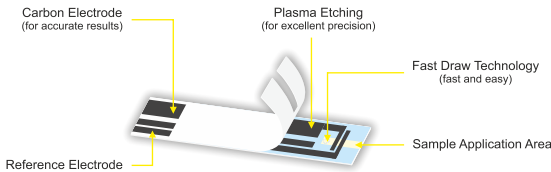

Specifications

Glucose ranges	:	20 - 800 mg/dL
Analysis time	:	5 sec
Sample volume	:	0.5 microliters
Sample	:	Whole blood
Sample mode	:	Capillary
Interference	:	Minimal
Fill detection	:	Automatic
Display	:	48 mm (Diagonally) LCD
Shelf life of strips	:	24 months (From Mfg. Date)
Strips packing	:	25, 50 and 100 (2 Vials x 50 Strips)
Memory	:	1000 Readings
Strips Storage Temp.	:	0-50°C
Operating Temp.	:	2-50°C
Relative Humidity	:	10-95% (Non Condensing)
Hematocrit	:	15-66%
Precision	:	+/- 4.5%
Battery	:	3v Li-ion Battery (CR2032), Lasts for 1000 tests

Soft Inserting Test Strip

Inserts the test strip easily and smoothly.

Gluco Spark is india's first, next generation, highly advanced ultra low volume glucose monitoring system available at affordable price. This glucometer uses highly advanced biosensor technology and requires 0.5 microliters sample volume and presents result in 5 sec.

Gluco Spark uses carbon electrode with to ensure precise & accurate results. **Gluco Spark** strips are manufactured with tight tolerances and do not require any coding.

Sensa Core's biosensor technology solves a number of immediate technical problems related to glucose strip by significantly enabling sensitive and accurate performance of glucose strip in a cost-effective way.